Mediblink[®]

medical devices

EN Mediblink Insect Bite Healer M002

INSTRUCTION MANUAL

Please read before use.

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Prije uporabe pročitajte.

TABLE OF CONTENTS

Intended use	7
Target group	7
Contraindication	7
Explanation of symbols	8
Safety instructions	10
General warnings	10
General precautions	11
Safety notes for batteries	12
Scope of supply	13
Power connection	14
Usage	14
What if there are problems?	17
Cleaning and care	18
Technical specifications	19
Product warranty	20
EMC information	22
Electromagnetic emissions	23
Electromagnetic immunity	24

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

EN Mediblink Insect Bite Healer M002

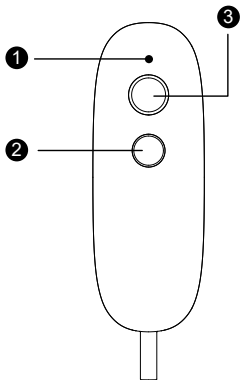- 1 Device status LED
- 2 Activation button
- 3 Power button

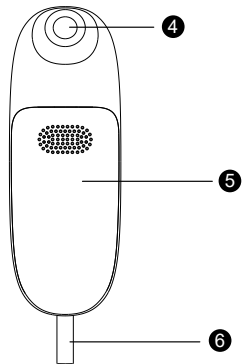- 4 Heating plate (MCH)
- 5 Battery compartment lid
- 6 Silicone rope

Intended use

The Mediblink Insect bite healer M002 is intended for relieving itching and swelling caused by insect bites or stings by applying heat (local hyperthermia).

Target group

The Mediblink Insect bite healer M002 is intended for use in the home healthcare environment. It can be used on anyone over 3 years of age under supervision, can be used for self-treatment by anyone over 12 years of age.

Contraindication

- Do not use the device if you suffer from diabetes. Sensitivity to pain could, in certain circumstances, be reduced and diabetics may burn themselves.
- Do not use the device on people with disabilities, children under 3 years old or people with reduced sensitivity to heat. This includes people with skin alterations due to illness or scarred tissue in the application area, after taking pain relief medication, drinking alcohol or taking drugs.
- Do not use the device on acutely or chronically diseased (injured or irritated) skin (e.g. inflamed skin, whether painful or not, reddened skin, rashes, (e.g. allergies), burns, frostbite, bruises, swellings (without signs of an insect bite/sting), both open and healing wounds, and surgical scars where the healing process could be affected).

- Do not use the device if you are suffering from a fever.
- Do not use the device on sensitive areas of skin.
- Do not use the device if you have any sensory impairment that reduces the feeling of pain (e.g. metabolic disorders).
- Do not use the device if you are using lotions, creams and gels locally at the same time.
- Do not use the device if you suffer from persistently irritated skin due to long-term heat application on the same area of skin.

Explanation of symbols

The following symbols are used on the device, in these instructions for use, on the packaging and on the type plate for the device:

	Caution
	Before use, consult instructions for use or electronic instructions for use: http://www.mediblink.com/f/m002.pdf
	Type BF Applied part
	Disposal in accordance with the Waste Electrical and Electronic Equipment EC Directive – WEEE
	Manufacturer

	Authorized Representative in the European community
	Date of Manufacture
	Serial Number
	Reference number
	Batch number*
IP42	Ingress protection
	CE Marking of Conformity, Notified body ID number
	RoHS mark
	Medical Device
	Unique Device Identifier

*The production date can be read from LOT number [YYYYMM]; first four digits represent the Year and last two digits Month of production. Example: LOT 202503 = March 2025

Safety instructions

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.

If you think you or someone in your family has experienced a serious reaction to a medical product, please report to manufacturer and the competent authority of the Member State in which you are established.

General warnings

- Adults and children aged 12 years and older may apply the product themselves; children under 12 years require application by an adult caregiver.
- Before using the device, consult your doctor if you have any skin conditions, serious illness (in particular if you have a propensity to thromboembolic conditions or recurrent malignant growths), unexplained chronic pain in any part of the body.
- Seek medical advice if symptoms persist despite using the insect bite healer, and discontinue use.
- Seek medical advice immediately if the symptoms worsen whilst or after using the device and stop using it.
- Seek medical advice immediately if it is a tick bite. The bite may result in transmission of pathogens such as tick-borne encephalitis (TBE) or borreliosis.
- Seek medical advice immediately if the first signs of increased itching, skin redness, swelling, dizziness, nausea, shortness of breath, formation of wheals or a drop or increase in blood pressure occur.

- Do not use near to the eyes, on the eyelids, mouth or on mucous membranes. Furthermore, the device must not be used on moist skin.
- Do not allow children to play with the device.
- Store the device out of the reach of children and pets.
- Keep packaging material away from children (risk of suffocation).
- Please note the general hygiene measures when using the device.
- The device is not a substitute for medical consultation and treatment. Always consult your doctor if you are experiencing any pain or are suffering from an illness.
- If the device does not work properly, or if you feel unwell or experience pain, stop using it immediately.

General precautions

- Do not touch the device with wet hands. Do not allow water to spray onto the device. Only operate the device if it is completely dry.
- Do not use any detergents or solvents to clean and maintain the device. If liquid penetrates the device this could cause damage to the electronics or other device parts and lead to a malfunction.
- The device should only be used by one person; use by multiple people is not recommended.
- Protect the device against direct sunlight and store in a dry place.
- Protect the device from impacts, moisture, dust, chemicals, marked temperature fluctuations and nearby sources of heat (ovens, heaters).
- Ensure that no water enters the device. If this should happen, remove the batteries and only use the device again when it is completely dry.

- If you store the device, no heavy objects should be placed on top of it.
- Before use, ensure that there is no visible damage to the device and karabiner. In case of any doubt, do not use the device, and contact your retailer or the specified Customer Service address.
- Store the device in a location protected against climatic influences. The device must be stored in the environmental conditions specified.
- When ambient temperature is less than $-10\text{ }^{\circ}\text{C}$ or higher than $50\text{ }^{\circ}\text{C}$, the time from minimum storage temperature to use is at least 4 hours; The time from maximum storage temperature is at least 4 hours.

Safety notes for batteries

- If your skin or eyes come into contact with battery fluid, rinse the affected areas with water and seek medical assistance.
- Choking hazard! Small children may swallow and choke on batteries. Therefore, store batteries out of the reach of small children.
- Risk of explosion! Do not throw batteries into a fire.
- If a battery has leaked, put on protective gloves and clean the battery compartment with a dry cloth.
- Do not disassemble, open or crush the batteries.
- Observe the plus (+) and minus (-) polarity signs.
- Protect the batteries from excessive heat.
- Do not charge or short-circuit batteries.
- If the device is not to be used for a relatively long period, take the batteries out of the battery compartment.

- Use identical or equivalent battery types only.
- Always replace all batteries at the same time.
- Do not use rechargeable batteries! Note on electromagnetic compatibility.
- The device is suitable for use in all environments listed in these instructions for use, including domestic environments.
- The use of the device may be limited in the presence of electromagnetic disturbances. This could result in issues such as error messages or the failure of the device.
- Avoid using this device directly next to other devices or stacked on top of other devices, as this could lead to faulty operation. If, however, it is necessary to use the device in the manner stated above, this device as well as the other devices must be monitored to ensure they are working properly.
- Failure to comply with the above can impair the performance of the device.

Scope of supply

Please check first of all that the unit is complete. The following parts are included as standard:

- 1 Insect bite healer
- 2 batteries (type AAA, LR03) 1.5 V
- 1 user manual

Please contact your supplier if you notice any transport damage on unpacking the unit.

Power connection

Inserting the batteries

1. Push open the battery compartment lid **5**.
2. Insert the batteries included in delivery into the battery compartment as shown inside the battery compartment. When inserting the batteries, ensure that the correct polarity (+/-) is observed.
3. Slide the battery compartment lid **5** closed until you feel and hear it click into place. Install the battery as above.

Remove the batteries if you do not intend to use the device for a long time (over 3 months).

Usage

WARNING

Before use, make sure that there is no insect sting remaining in the area that was stung. Remove the sting carefully before use. Use immediately after the insect bite or sting achieves the best results. If you wait too long to treat the area that has been stung or bitten and itching and swelling have already started, the insect bite healer can only have a limited effect. Nevertheless, treat the area with the insect bite healer. Usually, itching can be relieved and healing can be sped up.

1. Press the power button **3** to start the device. The device status LED **1** lights up green, a buzzer sounds and indicates to you that the device is ready to use.
2. For a 3-second application, place the heating plate (MCH) **4** of the device on the area that has been stung or bitten. Press the activation button once **2**. The device status LED **1** lights up blue and signals that the application is starting. The treatment temperature of the device is $47\text{ }^{\circ}\text{C} \pm 2\text{ }^{\circ}\text{C}$ ($116.6^{\circ}\text{F} \pm 3.6^{\circ}\text{F}$). This short application time is recommended for first use or use on sensitive skin.
3. For a 6-second application, place the heating plate (MCH) **4** of the device on the area that has been stung or bitten. Press the activation button **2** twice in quick succession. The device status LED **1** lights up yellow and signals that the application is starting. The treatment temperature of the device is $50\text{ }^{\circ}\text{C} \pm 2\text{ }^{\circ}\text{C}$ ($122^{\circ}\text{F} \pm 3.6^{\circ}\text{F}$). This application time is recommended for regular use on normal skin.
4. After the end of the application, a buzzer sounds, heating of the heating plate (MCH) **4** stops automatically and the device status LED **1** lights up green again. Remove the heating plate (MCH) **4** from the area of skin once the signal sounds.

IMPORTANT

Depending on the sensitivity of the person using the device and of the area of skin, the treatment temperature may feel unpleasant and cause slight skin redness and, in rare cases, irritated skin.

To achieve this itch-relieving effect on the skin nerves, a short period of heat pain is necessary. Only mild heat pain is required for treatment success. If the heat function feels too hot, stop treatment immediately.

5. If required, you can use the heat function of the insect bite healer on the sting/ bite again after a break of 2 minutes or also use it immediately on another sting/

bite in a different area. The maximum number of 5 applications per hour on the same treatment area must not be exceeded.

6. If the device is no longer required, press the power button **3** to turn off the device. The green device status LED **1** turns off. If there is no further operation within 5 minutes, the device will automatically shut down.
7. Before storing, make sure that the device is switched off and the heating plate **4** has cooled down.

Battery warning

Before the battery level becomes too low for other applications, the device status LED **1** starts to flash 3x and an acoustic signal sounds 3x. Replace both batteries soon.

What if there are problems?

Problem	Possible Cause	Solution
The device does not switch on.	The battery installation is incorrect.	Install the battery over again.
	The batteries are flat.	Replace both batteries as described in these instructions for use.
The device no longer gets hot enough.	The battery level is too low.	Replace both batteries as described in these instructions for use.
Device outputs 3 acoustic signals/the device status LED ❶ flashes – even after switching the device OFF and ON again.	Battery level is too low, meaning that the device may no longer be able to be operated.	Replace both batteries as described in these instructions for use.
Application cannot be started, even though the device status LED ❶ lights up green.	The device needs to be restarted.	Press the power button ❸ again.
Device outputs 5 acoustic signals/the device status LED ❶ flashes – even after switching the device OFF and ON again.	The device may be faulty.	Please contact Customer Services.

Cleaning and care

- The device is intended for multiple uses.
- It is recommended that the device is cleaned before first use and each additional use.
- The device must be switched off, allowed to cool down and the batteries removed each time before cleaning.
- Use only soft, dry cloths for cleaning the housing.
- Never use any substances for cleaning that could potentially be toxic if they came into contact with the skin or mucous membranes, or if they were swallowed or inhaled. Do not use any abrasive cleaning products and never submerge the device in water.
- Do not clean the device and karabiner in the dishwasher!
- Avoid mechanical cleaning with brushes or similar, as this could result in irreparable damage.
- When using the device on multiple areas that have been stung and bitten, it is recommended that you clean the heating plate each time. This reduces the possibility of germs and bacteria being transferred from one area of skin to the next. Use a commercially available alcohol disinfection wipe to disinfect the heating plate.

Disposal

This product must not be disposed of together with domestic waste.

All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner.

Remove the batteries before disposing of the device/unit. Do not dispose of old batteries with your household waste, but at a battery collection station at a recycling site or in a shop.

Consult your local authority or your supplier for information about disposal.

Technical specifications

Type	Insect bite healer
REF / Model	M002 / IBR10
Dimensions (L x W x H)	109.5mm x 34mm x 24.5mm
Weight (without batteries)	25g
Power supply	2 x batteries, type AAA (LR03 1.5V)
Max. temperature	50°C±2°C (122°F±3.6°F)
Functions	Sensitive application: 3s, 47°C±2°C (116.6°F±3.6°F) Regular application: 6s, 50°C±2°C (122°F±3.6°F)
Service life	5 years
Expected operating duration	Approx. 300 applications (per batteries)
Permissible storage and transport environment	-10°C to 50°C, relative humidity 30% to 85% pressure 700hPa to 1060hPa*
Permissible operating environment	10°C to 40°C, relative humidity 30% to 85% pressure 700hPa to 1060hPa

* In case of longer storage, remove the batteries and store the device in a cool place (below 30°C).

In accordance with our policy of continual product improvement, we reserve the right to make technical and visual changes without notice.

Product warranty

Product: Mediblink Insect Bite Healer M002

Manufactured for (Importer for EU and Distributor): Mediblink d.o.o., Gubčeva cesta 19, 8210 Trebnje, Slovenia; info@mediblink.com; www.mediblink.com

Seller's name, address, signature and stamp*:

Date of extradition/sales*:

*If the invoice is accompanied by this warranty, and if all above information can be seen from the invoice, it is not necessary to fill in this field.

Warranty terms

Dear customers!

The warranty period is **5 years** and starts from the day of product purchase. In case of product claim, you have to show the invoice. We kindly ask you to save the invoice!

Unfortunately, wrong handling with the device is a reason for 95% of customer complains. You can easily avoid any problem, by getting useful information provided by our special service department. To reach our service department, you can call or send e-mail to Mediblink local Distributor.

Before sending the product back to retailer, we kindly ask you to call our service department, to get help about how to use the device to save you with unneeded trips.

The Manufacturer guarantees free elimination of all imperfections due to defects in material or manufacturing procedure by repairing or replacing the product. In case that the product cannot be repaired or replaced, the customer will get the money refund. The guarantee is not valid in case of the force majeure, accidents or unexpected events (such as lightning, water, fire etc.), incorrect use or incorrect transport, non-compliance with safety and maintaining regulations or in case of unprofessional product intervention.

Traces of every day product usage (scratches, abrasions) and not subject to claim. The warranty does not eliminate the customer rights, which originate from seller responsibility for product flaws. By accepting the claimed product by the service department, the service department does not take responsibility for loss of saved data or settings on the product. All product repairs, which are performed out of product warranty period, have to be paid by customer by prior notice.

The Manufacturer guarantees the product quality and flawless product operation in the warranty period, which starts with the day of product purchase. If the product cannot be repaired in 45 days, the product will be replaced with a new one. In case that the product cannot be replaced, the money will be refunded to the customer.

In case of product claim, call or send e-mail to Mediblink local distributor.

EMC information

The ME EQUIPMENT or ME SYSTEM is suitable for home healthcare environments.

Warning: Don't use near active HF surgical equipment and the RF shielded room of an ME system for magnetic resonance imaging, where the intensity of EM disturbances is high.

Warning: Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.

Warning: Use of accessories, transducers and cables other than those specified or provided by the manufacturer of this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation.

Warning: Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the equipment, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.

The Mediblink Insect bite healer M002 / IBR10 may be interfered with by other equipment, even if that equipment complies emission requirements.

Electromagnetic emissions

Guidance and manufacturer's statement – Electromagnetic emission

Emission test	Compliance
RF emissions CISPR 11	Group 1
RF emissions CISPR 11	Class B
Harmonic emissions IEC 61000-3-2	Not applicable
Voltage fluctuations / flicker emissions IEC 61000-3-3	Not applicable

Electromagnetic immunity

Guidance and manufacturer's declaration – Electromagnetic immunity

Immunity test	IEC 60601-1-2 test level	Compliance level
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air
Electrical fast transient/burst IEC 61000-4-4	Not applicable	Not applicable
Surge IEC 61000-4-5	Not applicable	Not applicable
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	Not applicable	Not applicable
Power frequency magnetic field IEC 61000-4-8	30 A/m 50 Hz/60 Hz	30 A/m 50 Hz/60 Hz
Conducted RF IEC 61000-4-6	Not applicable	Not applicable
Radiated RF IEC 61000-4-3	10 V/m 80 MHz–2,7 GHz 80 % AM at 1 kHz	10 V/m 80 MHz–2,7 GHz 80 % AM at 1 kHz
NOTE U _T is the a.c. mains voltage prior to application of the test level.		

Guidance and manufacturer's declaration – Electromagnetic immunity

	Test Frequency (MHz)	Band (MHz)	Service	Modulation	Maximum Power (W)	Distance (m)	IEC60601-1-2 Test Level (V/m)	Compliance level (V/m)
Radiated RF IEC61000-4-3 (Test specifications for ENCLOSURE PORT IMMUNITY to RF wireless communication equipment)	385	380–390	TETRA 400	Pulse modulation 18 Hz	1.8	0.3	27	27
	400	430–470	GMRS 460 FRS 460	FM±5 kHz deviation 1 kHz sine	2	0.3	28	28
	710	704–787	LTE Band 13,17	Pulse modulation 217 Hz	0.2	0.3	9	9
	745							
	780							
	810	800–960	GSM 800/900, TETRA 800, iDEN 820, CDMA 850, LTE Band 5	Pulse modulation 18 Hz	2	0.3	28	28
	870							
	930							
	810	1700–1990	GSM 1800; CDMA 1900; GSM 1900; DECT; LTE Band 1, 3,4,25; UMTS	Pulse modulation 217 Hz	2	0.3	28	28
	870							
	930							
	2450	2400–2570	Bluetooth, WLAN, 802.11 b/g/n, RFID 2450, LTE Band 7	Pulse modulation 217 Hz	2	0.3	28	28
	810	5100–5800	WLAN 802.11a/n	Pulse modulation 217 Hz	0.2	0.3	9	9
	870							
930								

Guidance and manufacturer's declaration – Electromagnetic immunity

Radiated RF IEC61000-4-39 (Test specifications for ENCLOSURE PORT IMMUNITY to proximity magnetic fields)	Test Frequency (MHz)	Modulation	IEC60601-1-2 Test Level (A/m)	Compliance level (A/m)
	30 kHz	CW	8	8
	134.2 kHz	Pulse Pulse modulation 2.1 kHz	65	65
	13.56 kHz	Pulse modulation 50 kHz	7.5	7.5

Manufacturer:
 AOJ HEALTH Technology Co., Ltd.
 Building No. 1, 1-3F, Building No. 2, 2F,
 AOJ Industrial Park, Sanwei Community,
 Hangcheng Street, Bao' an District,
 Shenzhen 518126, China
 E-mail: marketing@aojmedical.com
 Made in China

**Manufactured for
 (Importer for EU & Distributor):**
 Mediblink d.o.o.
 Gubčeva cesta 19
 8210 Trebnje
 Slovenia
info@mediblink.com
www.mediblink.com

Share Info GmbH
 Am Schulzentrum 12, 41564 Kaarst,
 Germany
 E-mail: eu-rep@share-info.com

Instructions for use, version No: A0
 Issue date: 01.10.2023
 Date of last change: 30.11.2025

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